

July 2, 2021



Dear LaserCap Customer,

Your LaserCap kit includes a new version of the powerpack. Please refer to the following directions to charge your powerpack, and to operate your LaserCap.

Charging the Powerpack

1. Turn the powerpack switch to the charge (I) position.
2. Plug the wall charger into the black "IN" port of the powerpack.
3. The indicator light on the wall charger will be red when the powerpack is charging, and green when the powerpack is fully charged.
4. Disconnect the powerpack from the wall charger when fully charged.

Operating the LaserCap

1. Plug the yellow end of the power cable into the yellow "OUT" port of the powerpack.
2. Plug the black end of the power cable into the LaserCap device.
3. Place LaserCap on head.
4. Turn the powerpack switch to the discharge (II) position to begin treatment.
5. Note the green lights on the face of the powerpack indicate level of charge.
6. Turn the powerpack switch to the neutral (o) position to end treatment.


(OVER)

LaserCap Company
26 Alpha Park Drive
Highland Heights, OH 44143
1-855-424-7774
info@laser-cap.com
www.LaserCap.com

— THE ORIGINAL —
LASERCAP[®]




charge

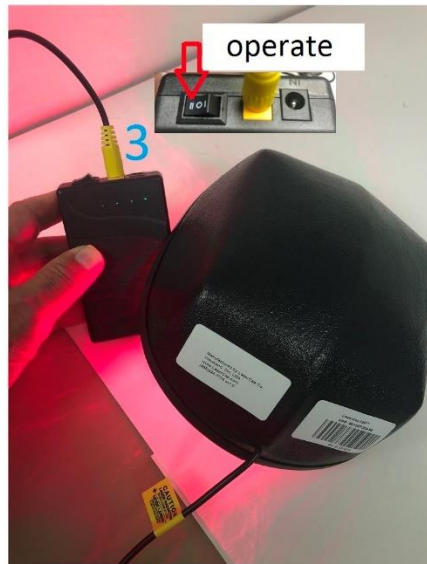

operate



charge-indicator

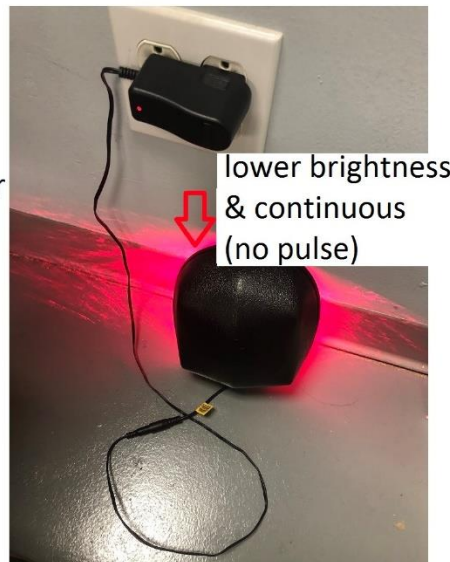

laser-cap-charger test

PowerPack (Model YB1202000) -

1. PowerPack 11.1V 2000mAh 22.2Wh
2. Charger 12.6V 1A Max 240V Input
3. Power Cable (Yellow End) 10Hz 50% Duty
4. BeltClip Pouch

Michael Rabin, MD
LaserCap Company
06-25-2021
rev4

LaserCap Company
26 Alpha Park Drive
Highland Heights, OH 44143
1-855-424-7774
info@laser-cap.com
www.LaserCap.com